



**Introducing The FORS™-15 Off-Loading Insole<sup>1</sup>**  
*The Durable and Adaptable Foot Off-Loading Solution*



**For Surgical Shoes**  
*(Shoe not included)*



**For Therapeutic Boots**  
*(Boot not included)*

**Consistent Off-Loading & Cushioning**

- Durable memory resistant Poron® for consistent cushioning, crush resistance, and long life.
- Highly stable off-loading zones; pre-cut plugs are torn out to create off-loading zones, but remain solidly in place when not disturbed. If you have ever experienced products with unstable and shifting pegs that produce inconsistent off-loading, you will welcome the durability and predictability of the FORS™-15 system!

**Best in Class Materials**

- Total of 15mm thick Poron® padding, with 5 mm always present under the foot, even after plugs are removed. Unique fabric mid-layer prevents off-loaded areas from collapsing. Tough and abrasion resistant Alcantara® top covers minimize slip and shear, absorb perspiration, are non-allergenic and bacteria resistant.
- Multilayered Poron® provides softness and support where needed: softer "slow rebound" Poron® adjacent to the foot, and denser Poron® below for support.

**Uniquely adaptable:**

- The FORS™-15 insole can be trimmed to any desired shape or size.

**Easy Care:**

- To completely clean and disinfect, soak the FORS™-15 insole in 10% chlorine bleach solution for 10 minutes and lightly scrub as needed, then thoroughly rinse and air dry.



3mm Extra-Soft Poron® Cushions foot


2mm Medium Density Poron® for Support

Fabric Layer Prevents Collapse and eliminates abrupt edges

>5mm Non-Removable Padding, Provides Constant Protection



**1** Alcantara® Cover, Absorbent and Shear Resistant

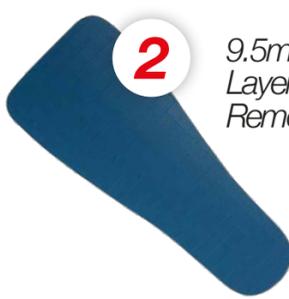

**2** 9.5mm Thick Poron® Layer, Pre-Cut with Removable 9.5mm Plugs

PRODUCT	DESCRIPTION	SIZES <sup>2</sup>	SHAPE
M494-15PZ	<i>FORS™-15 Insole: 15mm Multi-Layered Poron® Insole with Removable Plugs, Alcantara® Top Cover (Sized to fit Darco® MedSurg™ shoes<sup>4</sup>)</i>	<i>Women's Medium/Large<sup>3</sup> Men's Small Men's Medium Men's Large Men's Extra Large</i>	
M494-15PZ	<i>FORS™-15 Insole: 15mm Multi-Layered Poron® Insole with Removable Plugs, Alcantara® Top Cover</i>	<i>XXL (Fits XL Boots, can be trimmed as needed)</i>	

1. Patent Pending

2. FORS™-15 sizes are designed to correspond to industry standard surgical shoes and boots constructed with removable foam inserts. It is recommended to remove and replace the original foam insert with the FORS™-15 Insole. It is not recommended to place the FORS™-15 Insole on top of the original insole, as it may raise the heel out of the shoe or boot and create instability.

3. The Women's Medium size is created by trimming the heel of the Women's Large size. A trimming guide is printed on the insole.

4. Darco® and MedSurg™ are trademarks of DARCO International, Inc., 810 Memorial Blvd, Huntington, WV 25701, USA



**BONAPEDA Enterprises LLC**  
 Exclusive U.S. Distributor, Orthotics by Saluber

**www.bonapeda.com**  
**(973) 415-9090**